

说明书尺寸为:470x170mm, 图中比例为1:1

**料号:T-A050201 A版本**

**O'Yeet US**

Created in Pasadena, CA.  
Made for everybody.



@oyeetgo

# NEX Pro

USER GUIDE

**O'Yeet**

# BREAKTHROUGH YOUR LIMITS

We push the limits of technology to enhanced recovery and movement, anytime anywhere.

Powerful | Portable | Professional

## SAFETY FIRST

### SAVE THESE INSTRUCTIONS

We strive to make the NEX Pro as safe as possible. This is an advanced mechanical tool with electric components. If the NEX Pro and its accessories are not used or maintained properly, there is a risk of fire, electric shock, or injury. When using the NEX, these basic precautions should always be followed including:

#### USE ONLY AS INSTRUCTED.

Use the NEX Pro as described in this Instruction Manual. Use only NEX Pro recommended accessories and replacement parts. Do not carry out any maintenance other than as advised by NEX Pro.

#### NOT FOR CHILDREN.

The NEX Pro and charger are not intended for use by young children or persons with reduced physical, sensory or reasoning capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction by a responsible person. Do not allow the NEX Pro to be used as a toy. Close attention is necessary when used by or near children. Children should be supervised to ensure that they do not play with the NEX Pro or charger.

#### REDUCING RISK OF INJURY.

Charge the NEX Pro only in OYeet NEX Pro charging cable. Other charging cable may cause personal injury or damage. Do not connect the NEX Pro to a power supply plug or car cigarette lighter. The internal battery will be permanently disabled or damaged.

#### CHARGING LOCATIONS.

The NEX Pro should be charged indoors in a well ventilated, dry location. Do not charge the NEX Pro outdoors, in a bathroom or within 10 feet (3.1 meters) of a bathtub, or pool. Do not use the NEX Pro or charger on wet surfaces and do not expose the charger to moisture, rain, or snow. Do not use the NEX or charger in the presence of explosive atmospheres (gaseous fumes, dust, or flammable materials). Sparks may be generated, possibly causing a fire.

#### UNPLUGGING THE CHARGING CABLE.

Pull the plug, not the cord, to reduce the risk of damage to the electrical plug and cord. Keep the cord from heat, oil, and sharp edges. Do not use or charge with a damaged cord. Do not stretch the charger cord or place the cord under strain. Keep cord away from heated surfaces. Do not handle the Charging cable or NEX Pro with wet hands.

Unplug the charging cable when not in use. Unplug the charger from the outlet when not in use for extended periods. Make sure the charger cord is located so that it will not be stepped on, tripped over, or otherwise subjected to damage or stress. Have a damaged charging cable replaced immediately.

#### DO NOT BLOCK MOTOR VENTS.

Do not put any obstructions into vents. Do not use if any vent is blocked. Keep free of dust, lint, hair, and anything that may reduce airflow. Keep hair, loose clothing, and fingers from vents and moving parts.

#### DO NOT BURN OR INCINERATE THE NEX Pro.

The battery may explode, causing personal injury or damage.

Toxic fumes and materials are created when the battery is burned.

**DO NOT CRUSH, DROP, OR DAMAGE THE NEX CHARGING CABLE.**

Do not use a charger that has received a sharp blow, been dropped, run over, or damaged in any way.

BATTERY CHEMICALS CAUSE SERIOUS BURNS. Never allow the internal battery to come into contact with the skin, eyes, or mouth. If a damaged battery leaks chemicals, use rubber or neoprene gloves to dispose of it. If skin is exposed to battery fluids, wash with soap and water and rinse with vinegar. If eyes are exposed to battery chemicals, immediately flush with water for 20 minutes and seek medical attention. Remove and dispose of contaminated clothing.

**DO NOT OPERATE UNDER BLANKET AND PILLOW.**

Excessive heating can occur and cause fire, electric shock, or injury.

**STORING THE NEX Pro.**

Store in a cool, dry place. Only charge the NEX Pro when the ambient temperature is between 0°C/32°F and 40°C/104°F. Do not store the NEX Pro where temperatures may exceed 40°C/104°F such as in direct sunlight, in a vehicle, or metal building during the summer.

**DO NOT DISASSEMBLE.**

Disassembly or incorrect reassembly may result in the risk of electric shock, fire, or exposure to battery chemicals. The warranty will be void if the NEX Pro, or charger are disassembled or if any parts have been removed.

**SERVICE.**

If the NEX Pro is not working properly, has received a sharp blow, has been dropped, damaged, left outdoors, or dropped into

water, do not use and contact OYeet by emailing [help@oy-eet.com](mailto:help@oy-eet.com) for repair or service. Do not attempt to repair or disassemble the OYeet NEX Pro which may result in an electric shock or fire.

**DO NOT USE WHILE BATHING OR IN TUB OR SINK.**

Do not place the NEX Pro or charger where it can fall into a tub or sink. Do not place in or drop into water or other liquid. Do not reach for an appliance that has fallen into water. Unplug immediately.

**THERMAL LIMITER.**


The NEX Pro has an automatic reset thermal limiter that shuts off the NEX Pro to prevent overheating and fire.

# OYeet NEX Pro PRODUCT FEATURES


Model Name:	MG-S1
Frequency:	4 speed settings + AUTO rotate mode, 1600rpm-3500rpm,
Shock Depth:	12mm
Max Push Force:	>60lbs
Li-ion Battery :	16.8V    3*2550mAh
Charger Cable:	USB to DC 5V 1.5A
Continous Working Time:	>4 hours in normal load(below 5lbs) >1 hour in high load(5lbs - 12lbs)
Motor :	NEXPower™ QP38 by OYeet Dynamic Lab, 3 phases 14 poles with 3 sensors, rear earth magnets



# IN THE BOX

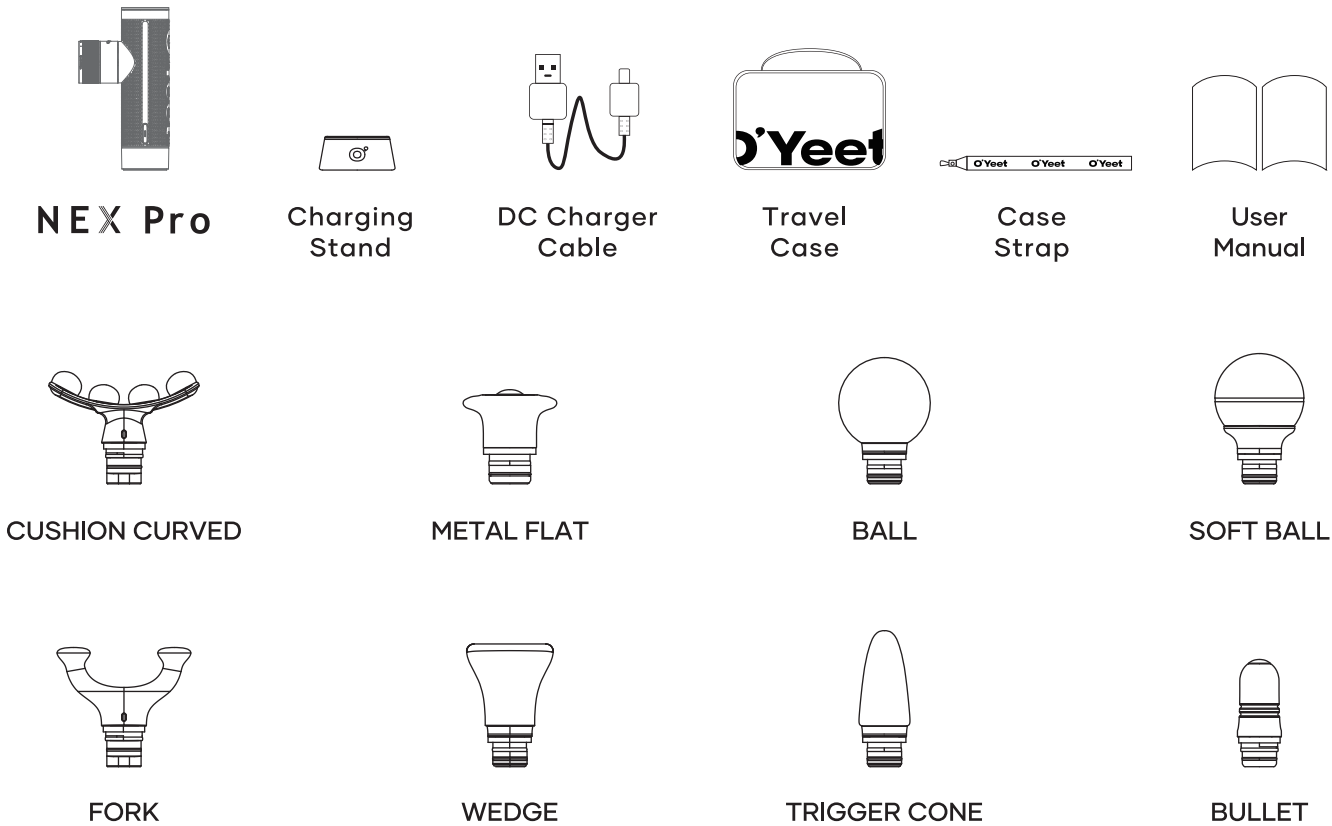


# KNOW YOUR O'Yeet NEX Pro

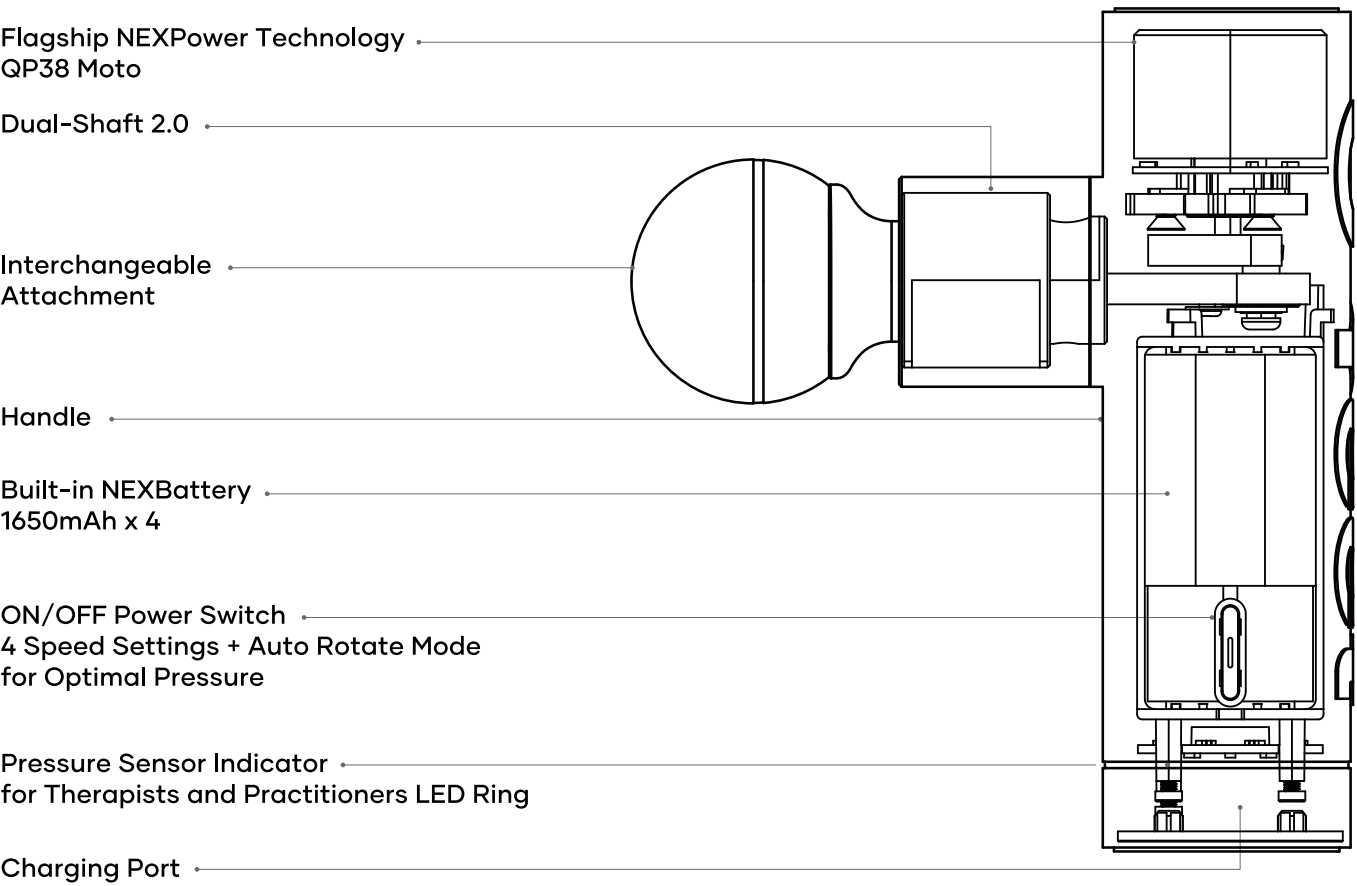

# USING YOUR OYeeet NEX Pro

- 1. To turn on the NEX Pro, firmly press and hold the ON/OFF power switch
- 2. Press the ON/OFF power switch to increase, decrease speed or set on Auto rotate mode.  
The NEX Pro has 4 speed settings and AUTO rotate mode which will rotate from 1600rpm to 3500rpm automatically.  
\* The LED light will start breathing when AUTO rotate mode is on.
- 3. To read your Pressure sensor indicator:  
The LED ring indicate your pressure level  
Yellow: 0 ~ 6.6lbs  
Orange: 6.6lbs ~ 13lbs  
Red: Greater than 13lbs
- 4. To turn off the NEX Pro, firmly press and hold the ON/OFF power button switch until the device shuts down.


**NOTE:** Fully charge NEX Pro for 6 hours before first use.

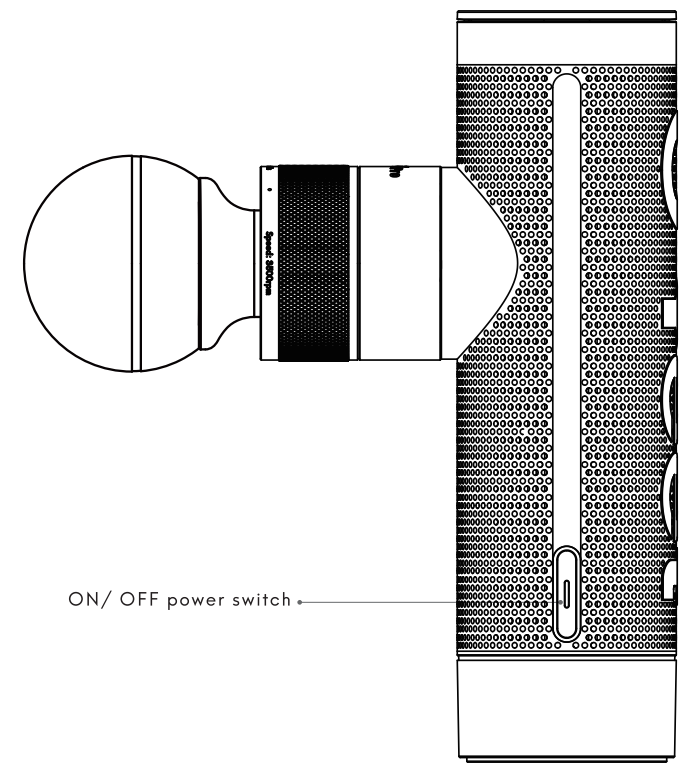

# CONNECTING ATTACHMENTS

- Align the attachment with the tip of the arm and push to connect.
- To remove, grip the attachment, placing 1 finger on each side of the plastic connector.
- Pull off the attachment.

 For attachment **FORK**, align the line on the attachment upward and push to connect.

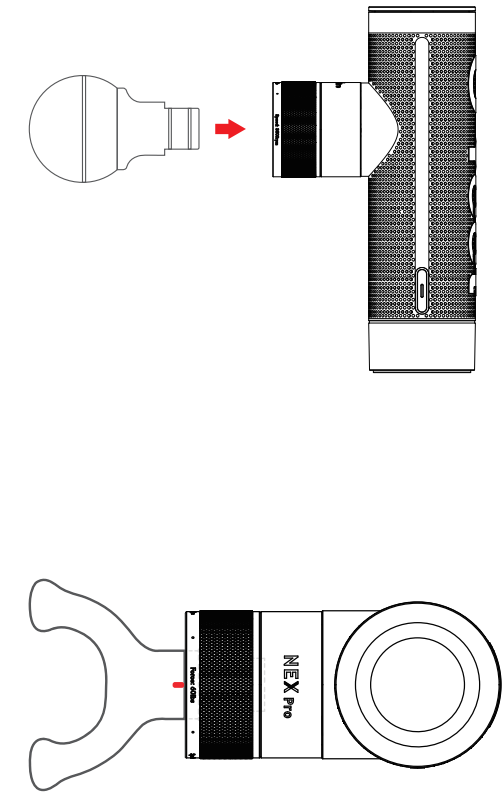

# CHARGING OYeeet NEX Pro

## CHARGING WITH DC CHARGER CABLE


The NEX Pro can be charged directly with DC Charger Cable. Simply connect the cable to the charging port and plug into power supply. The NEX Pro power switch LED light will turn off and the LED ring will start breathing when properly connected and stop breathing remain illuminated when fully charged.

## CHARGING WITH CHARGING STAND

Connect the DC Charger Cable to Wireless Charging Stand and plug into power supply, simply rest your NEX Pro on the stand and it will begin charging. The NEX Pro power switch LED light will turn off and the LED ring will start breathing when properly connected and stop breathing remain illuminated when fully charged.

## BATTERY

Green – 50% +  
Orange – 25~50%  
Red – 5~ 25%  
Flashing red light - below 5%, the device will be shuts down. Charge as soon as possible.

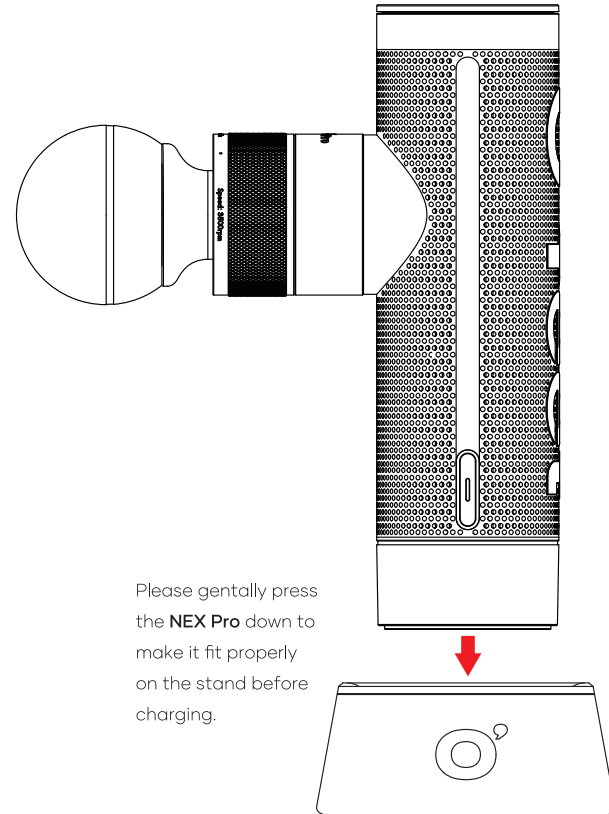

# ATTACHMENTS

The NEX Pro attachments set features 8 head attachments designed for specific treatments.


## METAL FLAT

Great to use on really dense muscle groups. Perfect for heavy practitioners.




## AIR CURVED

Simulated tapping massage simulating hand tapping massage with uniform intensity, such as arms, thighs and other curved parts.

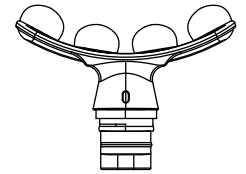


## WEDGE

Perfect for IT bands, shoulder blades and commonly used for "scraping" and "flushing", helping to decrease the amount of lactic acid in muscles.

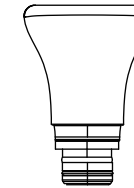


## FORK

Rounded ends to work around specific muscle groups.

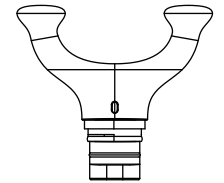

### SOFT BALL

Greater surface area, best attachment to start with to narrow down pain points.

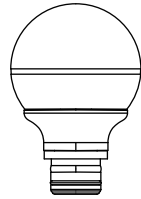


### BULLET

Pinpoint problem areas and release trigger points.



### BALL

Perfect for overall use to narrow down pain points.

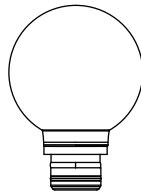

### TRIGGER CONE

Commonly used for pinpoint muscle treatment, including hands and feet.



## NEX Pro's ONE-YEAR EXTENDED WARRANTY

Simply register your product on **www.oyeet.com** within 14 days of your purchase to extend your warranty to **2 YEARS!**

## QUESTIONS?

Ask Us Anything!  
[help@oyeet.com](mailto:help@oyeet.com)